


Operating & Installation Guides

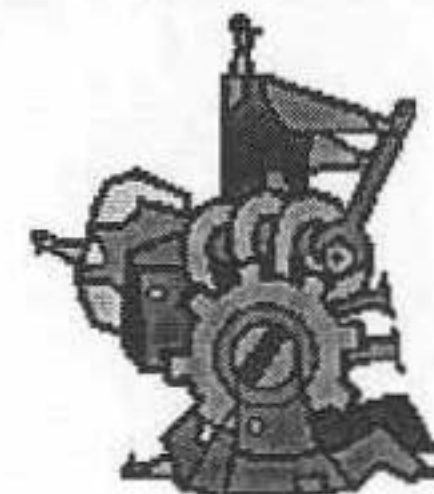

For the Vortex-Genie 2[®]

Operating and Installation instructions:

1. Plug the cord into a properly grounded, three-pronged outlet. (120V for the G560 model; 240V for the G560E; 100V for the G560-100V)
2. Your VORTEX-Genie[®] 2 is supplied with the 3" head installed and is ready for use for most one-handed applications. To install the cup head or other head attachments follow these steps:
 - a. Turn power to "OFF" (O) position
 - b. Remove the installed head by pulling straight up
 - c. Place alternate head on post and rotate until "flats" align
 - d. Press head down firmly until it snaps into place
3. For intermittent operation, push switch to "TOUCH" position. For continuous operation, push switch to "ON" position.
4. In either case, turn speed control to #3 in the shake range. (Mixer may not start if dial is lower than #3). In the "ON" position, you will see the head in motion immediately. In the "TOUCH" position, the vessel you are using must come in contact with the head to achieve motion. To get the desired mix, vary the speed of the motor and/or the angle.
5. When finished with either continuous or intermittent operation, return switch to "OFF" position in the center.
6. Your VORTEX-Genie[®] 2 mixer can become warm to the touch with constant use. There is nothing wrong with it. The stronger motor (UL listed) simply puts out more heat.

Caution! Do not lift your Vortex-Genie 2® mixer by its head.

All its heads, including the cup, are removable. They will pop off easily if you lift them, and snap on again just as easily by aligning and then pressing down.

Caution! Failure to follow operating instructions can compromise user's safety.

Care and Handling

Your VORTEX-Genie® 2 mixer should be given the care normally required for any electrical appliance. Avoid wetting or unnecessary exposure to fumes. The finish can be washed with water and soap or detergents, using a cloth or sponge.

Caution! Unplug from power before cleaning. Do not immerse.